

WarmTouch® Convective Air Warming Systems

CAREDRAPE® ADULT AND PEDIATRIC WARMING BLANKETS

CAREQUILT® ADULT AND PEDIATRIC WARMING BLANKETS

WARMTOUCH® 5200/5300 WARMING UNITS

Versatile blankets provide efficient heat delivery for optimal patient care.

CareDrape® Warming Blankets for the OR and CareQuilt® Warming Blankets for the PACU, achieve and maintain normothermia in your patients, which minimizes post-surgical complications and improves recovery time. The lighter, stronger material resists tears, punctures and fluids—providing more comfort and better care for your patients.

CareDrape Upper Body Warming Blanket



CareQuilt Adult Blanket (Pediatric Blanket available)



CareDrape Lower Body Warming Blanket



CareDrape Cardiac Blanket (Sterile)



Nellcor is committed to producing high-quality products that lead to better patient care.

Latex-free—No latex is used in any of our blanket fabrics, adhesives or packaging. No latex means no allergic skin reactions for your patient.

Unique design—Advanced quilted construction prevents billowing, allowing the blanket to conform to your patient's body.

Soft and comfortable—Fabric outer layers made of non-woven materials eliminate hot plastic surfaces that could contact patient's skin.

Economical—Convenient sizes to meet your clinical needs.

DESIGNED FOR COMFORT

CareDrape Warming Blankets

ALLOW EASY ACCESS TO PATIENTS.

CareDrape OR Warming Blankets for Upper and Lower Body and CareDrape Pediatric Warming Blanket

These blankets are specifically designed to withstand the demands of the OR environment. Made of a strong, two-ply material, they resist tears, punctures and fluids, yet are soft and flexible to permit custom draping. These three versatile blankets meet the needs of most surgical procedures, reducing inventory to save you money.


CareDrape Cardiac Blanket (Sterile)

This specialized blanket can be used for cardio-thoracic surgical procedures to regulate patient temperature peri-operatively. Two perforated slits allow easy access to the groin area, while the unique design minimizes air distribution at the patient's feet to maximize temperature regulation. A longer connector and inlet shield protects the sterile environment.

Rugged Construction

Fabric outer layer eliminates hot plastic surfaces that could contact patient's skin.

Layered construction adds strength, so even minor punctures won't expand, causing the blanket to deflate.

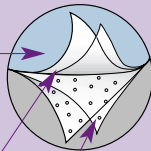

Precision die-cut holes help provide even airflow and a smooth, comfortable surface.

CareDrape Pediatric Warming Blanket

Position for under-body use

Immediate use in pre-op environment for patient comfort and warming prior to surgery.

Unique under-body design provides unrestricted access to surgical site.

Clear plastic cover traps warm air around patient; allows easy access and clear view of airway and IV lines.



Multiple inlets facilitate optimum blower placement.

CareDrape Upper Body OR Warming Blanket



X-ray transparent materials won't affect or interfere with X-rays.

CareDrape Lower Body OR Warming Blanket

Large size assures maximum skin surface coverage and warming.

Quilted design folds back easily for patient access.

Soft, flexible fabric allows custom draping for almost any procedure.



For more information call 1-800-NELLCOR,

WarmTouch 5200/5300 Warming Units For Use With CareDrape Warming Blankets


Fast warm-up, quiet operation and optimal temperature ranges for your OR or PACU needs. Assure patient safety and comfort when used with *CareDrape* and *CareQuilt* blankets.

Key Benefits

- Fast warm-up (38°C in 30 seconds)
- Choice of two models for PACU or OR use
- Multiple mounting options: Cart, Bedrail, IV Pole or Floor

5200 PACU Warming Unit

Temperature settings designed specifically for PACU use to maintain normothermia as patient recovers from surgery or post anesthesia.

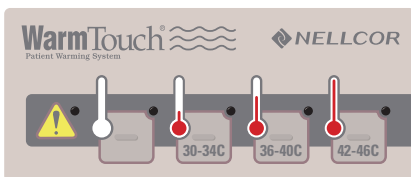

5200 PACU Control Panel

Automatic step-down feature from high setting (42°C-46°C) to medium setting (36°C-40°C) after 45 minutes promotes temperature management as the patient's thermoregulation stabilizes.


For more information or a demonstration contact your local Nellcor representative or call us at 1-800-NELLCOR (1-800-635-5267).



WarmTouch Warming Unit

5300 OR Warming Unit

Temperature settings designed specifically for use in the OR to offset patient hypothermia.

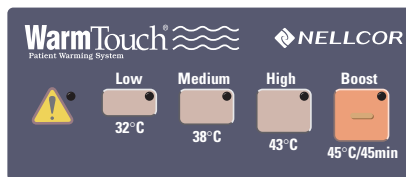

5300 OR Control Panel

- Ideal for the cold environments:
 - Preoperative area
 - Intra-operative (OR)
 - Post-operative
- Boost setting allows warming of patient at 45°C for 45 minutes
- Boost setting automatically reverts to high (43°C) after 45 minutes

The WarmTouch Warming Unit has multiple mounting options or may be placed directly on the floor.



Cart



Bedrail



IV Pole



Floor

Plastic warming cover (head drape) included.

Dual inlets facilitate optimum blower placement.

Medical-grade tape secures blanket and helps restrict airflow into surgical site.

WarmTouch Convective Air

WARMING SYSTEMS

Specifications

	CareDrape/Cardiac Warming Blankets	CareDrape Pediatric Warming Blanket	CareQuilt PACU Blankets
Size Uninflated:	Upper Body: 28 in W x 82 in L (71 cm W x 208 cm L) Lower Body: 41 in W x 56 in L (104 cm W x 142 cm L) Cardiac: 40 in W x 65 in L (102 cm W x 165 cm L)	25 in W x 41 in L (63 cm W x 104 cm L)	Adult: 40 in W x 79 in L (102 cm W x 165 cm L) Pediatric: 30 in W x 57 in L (76 cm W x 145 cm L)
Size Inflated:	Upper Body: 24 in W x 72 in L (61 cm W x 183 cm L) Lower Body: 34 in W x 52 in L (86 cm W x 132 cm L) Cardiac: 34 in W x 53 in L (86 cm W x 135 cm L)	22 in W x 35 in L (56 cm W x 89 cm L)	Adult: 33 in W x 74 in L (84 cm W x 188 cm L) Pediatric: 27 in W x 49 in L (69 cm W x 124 cm L)
Construction:	2-ply material consisting of polyethylene film inner layer and nonwoven outer layers.		
Material Characteristics:	Tear, puncture and fluid resistant; comfortable to the touch.		
Air Flow:	Quilted design uniform airflow distribution.		
Latex Allergen:	None.		

tyco

Healthcare

Nellcor

4280 Hacienda Drive
Pleasanton, CA 94588
Tel 925.463.4000
Toll Free 1.800.635.5267
www.nellcor.com

Tyco Healthcare UK LTD.
154 Fareham Road
Gosport, UK PO13 0AS
Tel +44.1329.224000

Specifications subject to change without notice.

© 2002 Nellcor Puritan Bennett Inc. All rights reserved.

B.I 0385-0602
WT06301

WarmTouch 5200/5300 Warming Unit

Size:	15 in H x 16 in W x 11 in D (38 cm x 41 cm x 28 cm)	AC Power:	115 VAC, 60 hz
Weight:	15 lb (7.3 kg)	Current Draw:	10 amperes, maximum
Hose:	5 ft (1.5 m); Collapsible and self-supporting	Heater Size:	1000 watts
AC Power		Warm-Up Time:	Air temperature at nozzle exceeds 38°C in approximately 30 seconds
Cord:	14-ft long (4.3 m); hospital grade	Air Flow Rate:	Approximately 55 CFM
UL Listing:	UL544, Hospital Medical Equipment	Noise Level:	Less than 55 db (normal conversation level)

Ordering Information

Description	Quantity	Catalog Number
CareDrape Upper Body Blanket with Head Drape	12/case	503-0870
CareDrape Lower Body Blanket	12/case	503-0880
CareDrape Pediatric Blanket	12/case	503-0850
CareDrape Cardiac Blanket (Sterile)	12/case	503-0860
CareQuilt PACU Blanket, Adult	25/case	503-0810
CareQuilt PACU Blanket, Pediatric	12/case	503-0840
WarmTouch 5200 PACU Warming Unit with Cart	1/box	501-5200
WarmTouch 5300 OR Warming Unit	1/box	501-5300
Replacement Filter for WarmTouch 5200/5300 Warming Unit	1/box	502-2200
5-Wheel Cart with Blanket Storage Basket	1/box	502-2900